


DATA SHEET

Deployment Manager

Manage tags and Exciters on-the-go from your iPhone®, iPad® or iPod touch®

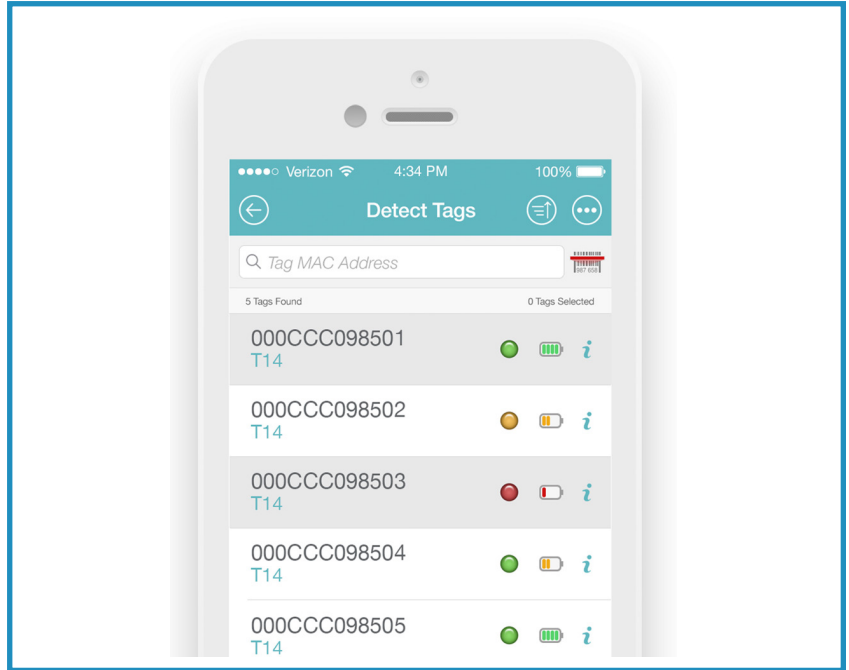

Product Highlights

- Connects via Bluetooth® to a TED device
- Connects and configures T15 sensor family tags using BLE technology
- Detect, activate and deactivate tags
- View tag details
- Configure bidirectional tags
- Reset tag battery counters
- Detect LF and US (Gen1 and Gen2 Exciters)
- View detected Exciter details
- Detect LF Interference

App Requirements

- TED Device: FW 1.3.5 or above
- Hardware Manager 2.1.20 or above
- iOS: 9.x, 10.x, 11.x and 12.x
- Supported Devices:
 - » iPhone 5S and above
 - » iPad Pro, Mini, Air
 - » iPad 3rd generation and above
 - » iPod touch 5th generation and above

Description

Deployment Manager is a component of STANLEY Healthcare’s visibility platform, which connects to TED devices and T15 family tags via Bluetooth® Low Energy (BLE) technology. The app can be used to detect, activate and view tag details, configure bidirectional tags, detect and view Exciter details, and detect LF interference.

The app is designed for use on the iPhone®, iPad® or iPod touch® and can be downloaded from the App Store®.



Key Features

Tag Management

Tag Management allows you to perform the following:

- Detect, activate and deactivate tags
- View tag details
- Configure bidirectional tags (T12, T14, P10 and T15 tag families)
- Configure sensor tag settings (T15 sensor tag family)
- Initialize tag batteries (reset battery counters)
- Make tags blink (supported tags only)

Detect Exciters

Detect Exciters allows you to perform the following:






- Detect LF and US (Gen 1 and Gen 2) Exciters
 - » Shows good LF and US signals
 - » Shows bad LF signals
- View Exciter details
- View Exciter RSSI history graph

Interference

- Detect LF interference

TED Device

- View current status and details of the connected TED device

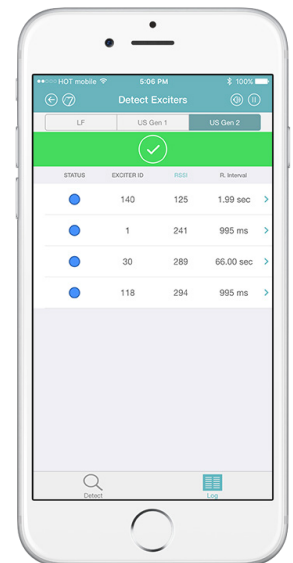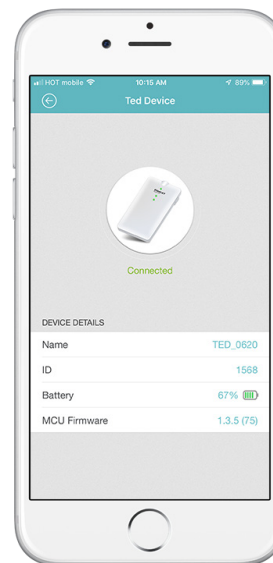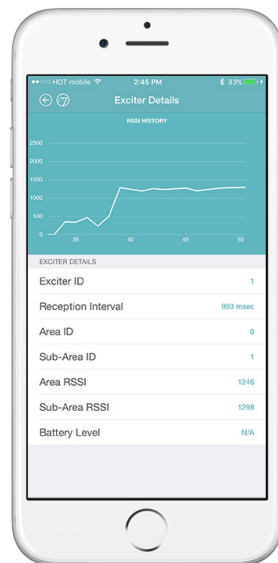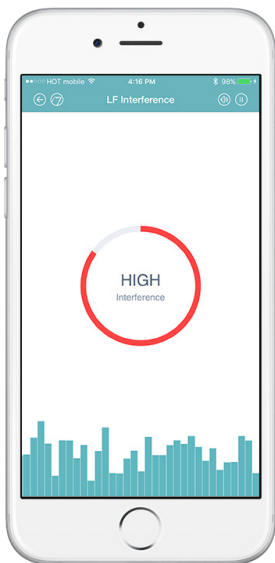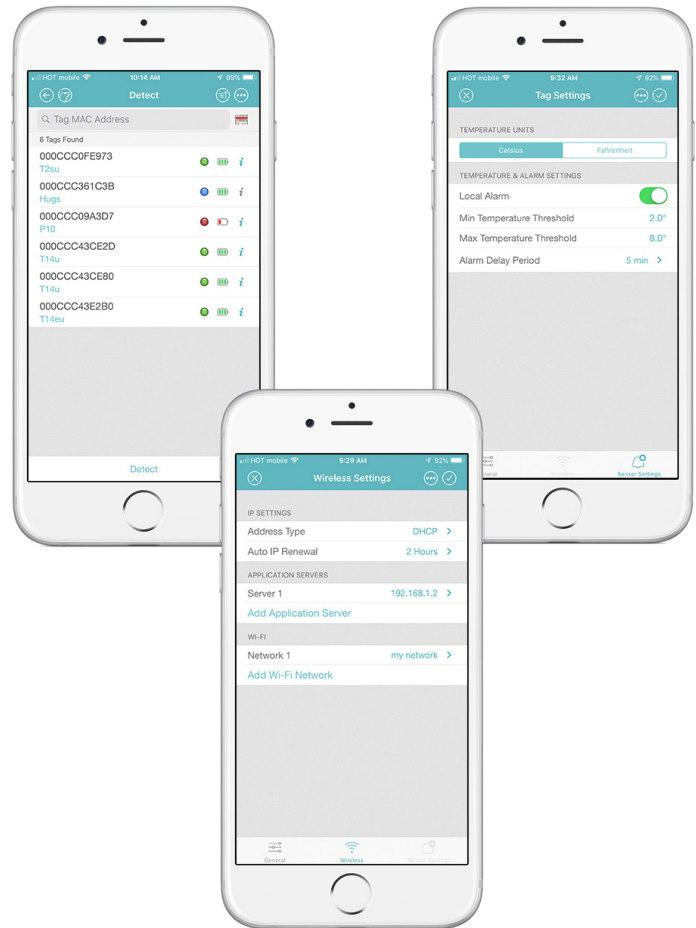

About STANLEY Healthcare

Over 15,000 hospitals and senior living communities rely on STANLEY Healthcare solutions to empower caregivers to deliver better care. STANLEY Healthcare is a part of Stanley Black & Decker and a proud supporter of the Alzheimer's Association®. Learn more at www.stanleyhealthcare.com.

130 Turner Street, Suite 700, Waltham, MA 02453 | 1-888-622-6992 | stanleyhealthcare.com | stanleyhealthcare@sbdinc.com